

Vacuum Filter Funnel System

Description

This Vacuum Filter Funnel is made of durable polypropylene and stainless steel. It is easy to set up and use and consists of:

Vacuum Filter Funnel

Description	Qty	Part No.
Alltech® 1L Vacuum Filtration System	Ea	6409
Alltech® 2L Vacuum Filtration System	Ea	6429

Replacement Parts for 1L System

Description	Qty	Part No.
1L Replacement Funnel	Ea	6473
Hose Barb Adapter Kit (with 8 adapters)	Ea	6438
Stainless Steel Mesh Filter Support Ring	Ea	6421
Silicone Rubber O-Ring	Ea	8619230
Retainer Ring	Ea	8619231

Assembly

Assemble the unit for use as follows:

1. Secure Base Unit (E) to receiving bottle. The receiving bottle must have a GL-45 threaded neck. The receiving bottle should never exceed 2 liters in volume. It is recommended that safety coated receiving bottles be used.
2. Install the Stainless Steel Membrane Support (D) into the Base Unit (E). This will fit snugly.
3. Apply the 47mm Membrane Filter (C) to rest on the Stainless Steel Membrane Support.
4. Install the Filter Retainer Ring (B). This ring protects the Membrane Filter and enhances the seal between the Funnel and Base Unit. The unit may be used without the Filter Retainer Ring, but great care must be taken to avoid damaging the Membrane Filter.
5. Seat the Funnel (A) on top of the Filter Retainer Ring (B) and raise the white cap to engage the threads on the Funnel. It is recommended that the cap be turned to tighten and seal the Funnel to the Base Unit, rather than turning the Funnel.
6. Secure the unit with a ring stand. A three-finger clamp can be used to hold the base unit, securely fixing the entire filtration apparatus to the ring stand.
7. Thread the 1/4"-28 to female Luer adapter into the vacuum port. Apply the appropriate barb to interface the unit to the vacuum source.
8. Apply vacuum and use the filter.

Safety Precautions

1. Never pull vacuum on any container unless it is designed and rated for vacuum service.
2. Protect against failure and do not use equipment beyond its design limits.
3. Ensure fluid compatibility with the filter system elements.
4. Always wear proper personal protective equipment appropriate to the task performed.

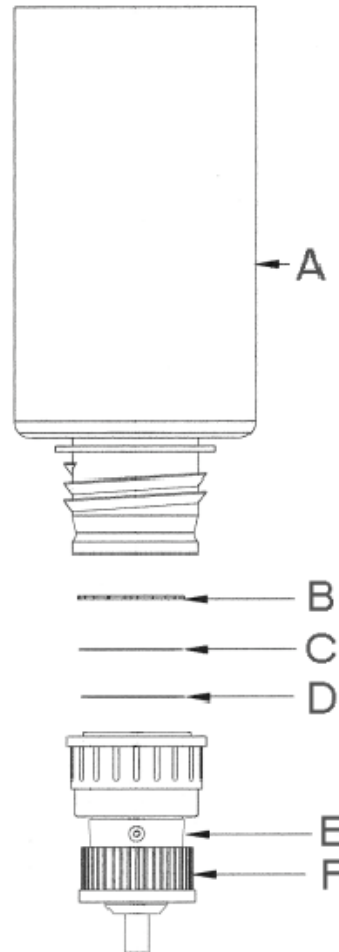


Figure 1 – Components of the Vacuum Filter Funnel and their proper orientation and assembly.

- A – Funnel
- B – Filter Retainer Ring
- C – 47mm Filter Membrane
- D – Stainless Steel Membrane Support
- E – Base Unit with Vacuum Port
- F – GL-45 Cap Closure

www.discoverysciences.com

Grace Davison Discovery Sciences • 2051 Waukegan Road • Deerfield, IL 60015
 Telephone: 847.948.8600 • Fax: 847.948.1078 • E-mail: DiscoverySciences@grace.com

2/2010

ALLTECH® is a trademark, registered in the United States and/or other countries, of Alltech Associates, Inc. GRACE® and GRACE DAVISON® are trademarks, registered in the United States and/or other countries, of W. R. Grace & Co.-Conn. GRACE DAVISON DISCOVERY SCIENCES™ is a trademark of W. R. Grace & Co.-Conn. This trademark list has been compiled using available published information as of the publication date of this brochure and may not accurately reflect current trademark ownership.

Alltech Associates, Inc. is a wholly-owned subsidiary of W. R. Grace & Co.-Conn. Grace Davison Discovery Sciences is a product group of W. R. Grace & Co.-Conn., which now includes all product lines formerly sold under the Alltech brand. © Copyright 2010 Alltech Associates, Inc. All rights reserved.

The information presented herein is derived from our testing and experience. It is offered for your consideration and verification. Since operating conditions vary significantly, and are not under our control, we disclaim all warranties on the results that may be obtained from the use of our products. Grace reserves the right to change prices and/or specifications without prior notification.